


Bari Foam+



Width adjustable from 915mm (36") - 1220mm (48") using foam inserts

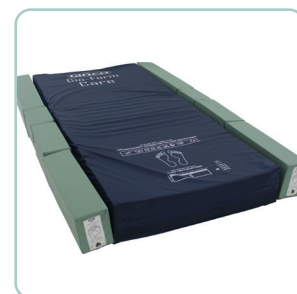

Enclosed in PU coated box cover

The Bari Foam+ mattress is an alternative to our dynamic air mattress and provides pressure relief for medium risk patients.

The mattress is unique with its width adjustable feature to compliment the aurum bariatric bed. Width adjustment is achieved through two foam extension inserts which slot either side of the mattress when the bed is extended. These foam extension pieces also provide a more solid edge for patient ingress and egress.

The mattress is suitable for patients up to 320kg (50st) and is made up of a series of combustion modified foam to give the mattress its pressure relieving qualities.

Features	Benefits
320kg (50st) swl	Able to cater for a range of patient weights
Width expandable from 915mm (36") - 1220mm (48") using two foam extension pieces	Compliments the width adjustment of the aurum bariatric bed
PU coated top cover	Vapour permeable and waterproof
Made up of a series of combustion modified foam layers	Provides pressure relief for medium risk patients
Removable washable cover	Infection control compliant

Product codes

BMPB065FP Bari Foam+ Bariatric Foam Mattress